


V-Block™

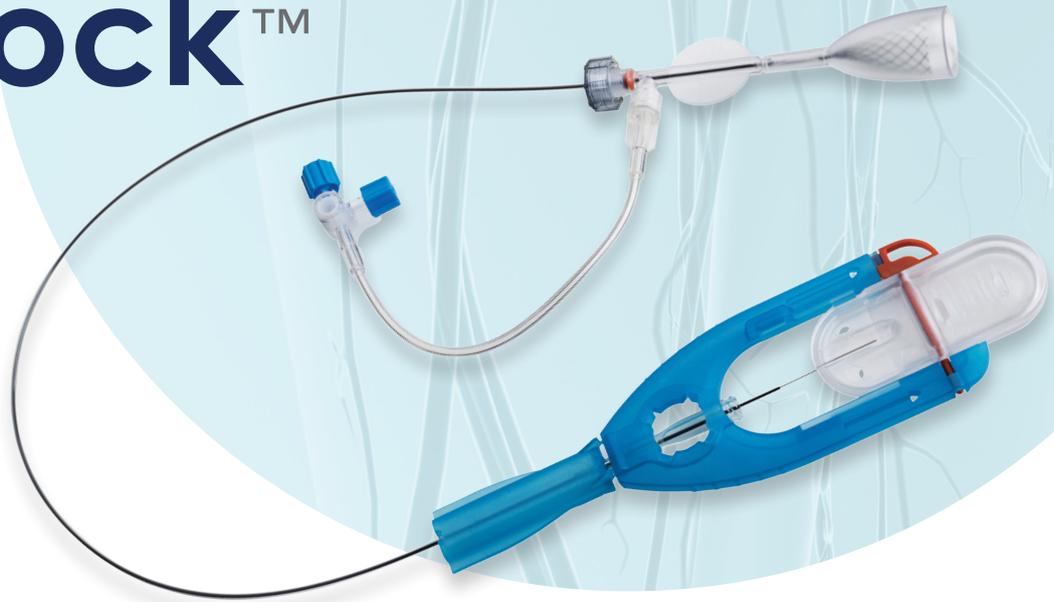

ScleroSafe with V-Block is especially designed to simply, safely and effectively treat severely dilated and challenging Great Saphenous Vein (GSV) conditions, including GSV valve incompetence and Reflux.

While all ScleroSafe products share the Double Lumen Catheter and Inverse Action Dual Syringe design concepts, V-Block also features a unique embolization device which further enhances the safety and efficacy of treating larger vessels.

Indicated for the occlusion of lower extremity veins with a diameter in the range of 4 to 14mm.

ScleroSafe with V-Block can be used with the 10cc/20cc Dual Syringe, a proprietary, specially designed (6F/450mm or 6F/650mm) catheter and a dedicated guide wire.

The embolization safety device is designed to block the blood flow within GSVs of larger diameters during ECA/Sclerotherapy procedures. The purpose of the device is to minimize the risk of upflow of the injected substance to the femoral vein and to the upper body. It is a biomedical-grade, flexible nitinol frame covered with a polymer-coated mesh and is tested and approved as a permanent venous implant.



Order Code	Description	Master Pack
VBL0001	V-Block 45 cm	10 units
VBL0002	V-Block 65 cm	10 units

Inverse Action Dual Syringe

(10ml / 20ml)



The unique inverse action dual syringe (also known as the DPS) injects the ECA (Endovenous chemical ablation) substance into the vein through one syringe while simultaneously aspirating blood through the second syringe – all within a single motion and with single-hand operation.

The aspiration of the blood results in negative pressure just behind the injected sclerosant agent. This key feature prevents uncontrolled overspill of the sclerosant upstream, ensures that there is no dilution of sclerosant with blood and minimizes the risk of inflammation. There is minimal pain, no scarring, discoloration or matting, no trapped blood clots, minimal to zero thrombophlebitis and immediate, visible results. VVT's unique syringe is compatible with both foam and liquid ECA substances.

For use with the V-Block.

Order Code	Description	Master Pack
DPS-2010	Inverse Action Dual Syringe (DPS) 20ml/10ml	32 units